

HeartSave myPAD

IT'S IN YOUR HANDS



Quick user guide

Turning on AED and attaching electrodes

1




To switch on the AED, press the green button until the green LED lights up and you hear voice instructions.

2




Call emergency services!

3

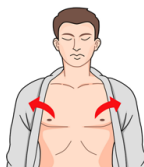


Remove clothing from the patient's upper body.
Cut clothing if necessary.

4

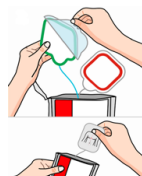

Open the electrode package, remove the protective foil.
Use the feedback sensor if available.

5

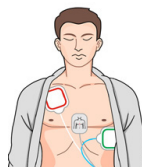

Attach the electrodes to the patient's bare upper body as shown.
If necessary, plug the connector into the AED.

Semi-automatic mode

Full automatic mode

6



Do not touch the patient!
The AED will automatically start the ECG analysis.



7




If defibrillation is required, the shock button will light up and a loud warning sound will be heard.




If defibrillation is required, the shock indicator will light up.



8

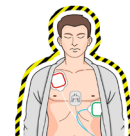

Make sure that neither helpers nor bystanders are touching the patient!

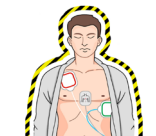

9



Press the shock button to release the shock.



The AED will count down a 3-second timer and then automatically initiate defibrillation.



10




Follow the AED's instructions for CPR.



CPR


Laypersons

11

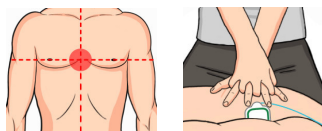

2 minutes chest compressions in time with the metronome.
100 per minute, 5-6cm deep.

Trained first responders

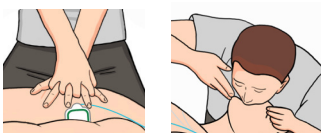

2 minutes alternating:

- 30 chest compressions
- 2 rescue breaths

12



Afterward, the AED will start another ECG analysis.

Do not touch the patient during this!

13



Continue all measures until professional help arrives.